MEYRA®

iMOVE

stay active!

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iMOVE

stay active!

THE ACTIVE FOOTREST FOR MORE MOVEMENT WHILE SEATED

EFFECTIVE. CAN BOOST CIRCULATION AND METABOLISM, LOOSENS AND RELAXES

EFFORTLESS. MOBILISES VIRTUALLY THE ENTIRE BODY WITHOUT ANY EFFORT AS A SIDE EFFECT

FLEXIBLE. COMPACT AND EASY TO TRANSPORT, QUICKLY READY FOR USE

SIMPLE. INTUITIVE AND VERY EASY TO USE

ENGLISH

CONTENT

1 USER SAFETY	18
2 SAFETY PROVISIONS	19
2.1 Hazards posed by electrical power	19
2.2 Structural modifications	19
2.3 Copyright	19
3 SCOPE OF DELIVERY	19
4 INITIAL SETUP	20
4.1 Positioning and connection	20
4.2 Operating modes	20
4.2.1 Operation directly on the device	20
4.2.2 Operation using the app	21
4.2.3 Switching between "smartphone" and "remote control" mode	22
4.2.4 Installing the remote control	22
4.2.5 Remote-controlled operation	22
5 RECOMMENDATIONS FOR USE	23
5.1 First steps	23
5.2 Suggested exercises	23 - 24
6 SPECIFICATIONS	24
6.1 Dimensions and weight	25
6.2 Electrical data	25
6.3 Marking	25
6.4 Material	25
7 CARING FOR YOUR iMOVE	26
8 GUARANTEE POLICY	26
9 PROBLEMS AND HOW TO FIX THEM	27
10 MAINTENANCE AND CUSTOMER SERVICE	27
11 DISPOSING OF YOUR iMOVE	27

1 USER SAFETY

Please read the user manual through carefully before you start using the device. The user manual contains important information on how to use your iMOVE. MEYRA GmbH shall not be liable in the event of the user failing to observe instructions and warnings. Keep this user manual in a safe place for future reference. Check the device for damage on receipt. In the event of damage, do not start using the iMOVE! If in doubt, please contact MEYRA GmbH.

POWER CONNECTION

Check whether the values indicated on the nameplate of the iMOVE correspond with the power supply of the location where the device is being placed.

GENERAL

- Please read the user manual thoroughly before using the iMOVE.
- Should any problems arise while using the device, please switch it off immediately and contact the MEYRA GmbH customer service centre.
- The iMOVE is not intended for use in potentially explosive locations.
- The iMOVE is not water-resistant and must not come into contact with water.
- In order to ensure ventilation, the iMOVE must not be covered whilst in operation.
- The iMOVE is for indoor use only.

GUARANTEE AND LIABILITY

Our "General terms and conditions" shall apply in all cases. In the event of physical injury or damage to property, no claims in respect of guarantee or liability can be accepted if such injury or damage can be attributed to one or more of the following causes:

- misuse
- opening of the device or improper initial setup and operation
- unauthorised structural modifications
- improperly executed repairs
- disasters due to external influences and force majeure

2 SAFETY PROVISIONS

2.1 HAZARDS POSED BY ELECTRICAL POWER

No hazards are posed by electrical power when the device is in perfect technical condition.

DANGER

Should the cable be faulty, it must be replaced immediately to avoid the risk of an electric shock during contact.

2.2 STRUCTURAL MODIFICATIONS

No modifications, additions or alterations may be made to the device without the prior consent of the manufacturer. Any such measures require confirmation in writing from MEYRA GmbH. Repairs may only be carried out by companies licensed by MEYRA GmbH.

2.3 COPYRIGHT

MotoSana GmbH retains the copyright to this user manual. This user manual is intended for the user only. It contains provisions and notes which may not be reproduced, distributed or otherwise communicated either in full or in part.

3 SCOPE OF DELIVERY

- Base unit iMOVE
- Power supply unit and connector plug
- Remote control (optional)
- User manual

4

INITIAL SETUP

4.1 POSITIONING AND CONNECTION

A level and low-vibration surface is important to ensure the iMOVE operates smoothly. The iMOVE is only authorised for indoor use. Having properly positioned the iMOVE, insert the connector plug into the socket. Once positioned and connected, the device is immediately ready for use.

CAUTION!

The device must not be used in wet rooms. Temperature range from 0 to +30 degrees Celsius. The device switches off automatically around four hours after the last command.

4.2 OPERATING MODES

The iMOVE can be operated directly on the device, by remote control and from your smartphone or tablet with the free iMOVE app.

4.2.1 OPERATION DIRECTLY ON THE DEVICE

There are three buttons on the left-hand side of the iMOVE with which you can operate the device.

1 STARTING AND STOPPING

Press to start and stop the iMOVE.

When is pressed, the iMOVE automatically starts with a frequency of 2 hertz.

Pressing again switches the iMOVE off. The backlights of the and buttons go out.

**2 SETTING THE FREQUENCY**

Press to increase the frequency during the exercise programme.

Press to decrease the frequency during the exercise programme.

4.2.2 OPERATION USING THE APP

The iMOVE app is available for iOS and Android. Download the app free of charge on to your smartphone or tablet from App Store (iOS) or Google Play Store (Android). You first have to pair your smartphone or tablet with the iMOVE via Bluetooth. In order to do so, check which mode the iMOVE is in. If the button is lit up blue, the iMOVE is in "smartphone" mode. If the button is lit up white, the iMOVE is in "remote control" mode. (To switch between "smartphone" and "remote control" mode see 4.2.3.) The button lighting up blue indicates that the iMOVE is ready to pair with your smartphone. Start the iMOVE app on your smartphone or tablet and open the Bluetooth settings at the top right in the iMOVE app. Select your iMOVE to pair the device.

1 STARTING AND STOPPING

Pressing the START button automatically starts the exercise programme with a frequency of 2 hertz.

Press PAUSE to halt the exercise during use and continue it at a later time.

Pressing STOP ends the exercise.

2 SETTING THE FREQUENCY

You can adjust the frequency either prior to starting the exercise programme or during the exercise.


Use the slider or the plus and minus buttons to increase or decrease the frequency.

3 SETTING THE STOCHASTIC

Use the stochastic slide button to switch the stochastic on and off. Stochastic means that the frequency changes at random in a range from -2 to +2 units in relation to the selected frequency.

4 EVALUATION

Having completed the exercise, a brief evaluation of your exercise unit is displayed.

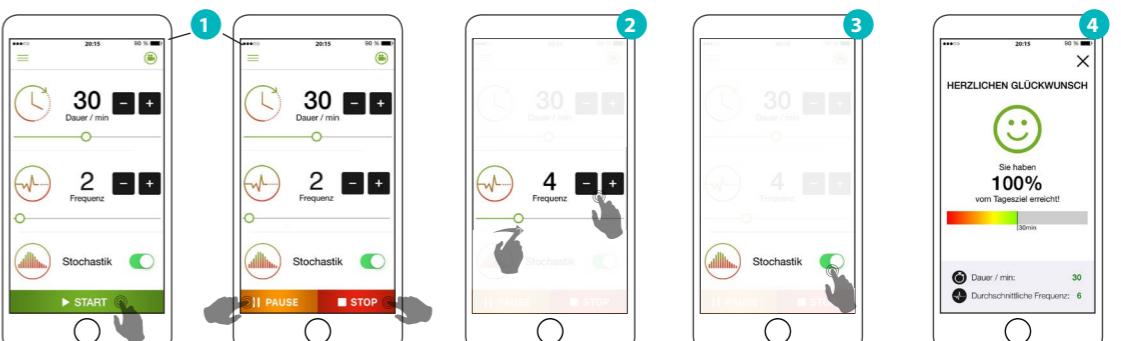4.2.3 SWITCHING BETWEEN "SMARTPHONE" AND "REMOTE CONTROL" MODE

The colour of the LED backlight of the button indicates the current mode of the iMOVE. If the button is lit up white, the iMOVE is in "remote control" mode. If the button is lit up blue, the iMOVE is in "smartphone" mode.

To switch between "remote control" and "smartphone" mode, proceed as follows:

Switch the iMOVE off. Press the button and the button simultaneously and keep both buttons pressed for approximately three seconds. The backlight of the button changes colour and the iMOVE is now in the corresponding mode. You may switch as required between the two operating modes.

4.2.4 INSTALLING THE REMOTE CONTROL (OPTION AVAILABLE AT A SURCHARGE)

If you ordered your iMOVE with remote control (optional), this is already paired with the device. To subsequently pair the remote control with the iMOVE, proceed as follows:

- 1 Check which mode the device is in. The iMOVE has to be in "remote control" mode.
See point 4.2.3
- 2 Press the button on the device and keep it pressed for around five seconds until it starts blinking blue.
- 3 Now press the button on the remote control for around five seconds until a green light lights up briefly.
The remote control is now connected with your device.

4.2.5 REMOTE-CONTROLLED OPERATION (OPTION AVAILABLE AT A SURCHARGE)

1 STARTING AND STOPPING

- Press to start and stop the iMOVE.
- When is pressed, the iMOVE automatically starts with a frequency of 2 hertz.
- Pressing again switches the iMOVE off.

2 SETTING THE FREQUENCY

- Press to increase the frequency during the exercise programme.
- Press to decrease the frequency during the exercise programme.

You can pair your iMOVE with four devices in total (smartphone, tablet, remote control). In the event of a fifth pairing, the oldest pairing to date is deleted. Press the button and the button simultaneously and keep both pressed for approximately ten seconds. When all three buttons blink briefly, this indicates that all pairings are deleted.

5

RECOMMENDATIONS FOR USE

5.1 FIRST STEPS

Start with low frequencies. Sit on a chair and place your feet on the pedals of the iMOVE. Switch the iMOVE on. Vary the speed to get a feel for the iMOVE. Find the intensity that feels best for you. Try out various foot positions. The exercising effect varies depending on the foot position. Also vary the counter-pressure you are exerting and change the position of your body. By resting your lower arms on your knees you can exercise your shoulders and the upper part of your body at the same time.

5.2 SUGGESTED EXERCISES

The iMOVE is flexible, compact and quiet - the perfect footrest at home, in the office or out and about while on holiday or travelling. It enables uncomplicated exercising and is suitable for virtually everyone. Reports based on experience using the iMOVE show: Exercising with the iMOVE activates the muscles, stimulates metabolism, improves circulation, and loosens and relaxes the muscular system. The iMOVE gets you moving while protecting your joints. Exercising with the iMOVE benefits the elderly and physically impaired as well as children and athletes alike. Use it while occupied with other activities like reading, watching television or working.

IMPORTANT!

Should any physical complaints arise, please stop exercising and consult your doctor.

**TO SOOTHE TIRED AND HEAVY LEGS**

Place both feet on the iMOVE and lean back. Select the speed that suits you best and relax.

**SHEER RELAXATION – TO HELP YOU GET TO SLEEP**

Lie down in front of the iMOVE and lightly rest your heels on it. Select a slow speed.

**TO LOOSEN SHOULDERS AND NECK**

Place your feet slightly further forward on the iMOVE. Gently rest your elbows on your knees to direct the vibrations into your shoulders and back.

**TO LOOSEN SHOULDERS AND NECK**

Position yourself on all fours in front of the iMOVE. Place your hands on the iMOVE (one hand on each pedal) and select a medium speed.

**TO IMPROVE CIRCULATION**

Place both feet on the iMOVE and roll from your heel to the ball of your foot in a seesaw movement. Select a moderate speed.

**TO IMPROVE CIRCULATION**

Place both feet on the iMOVE and roll your feet from the outer to the inner side of the foot and back.

6**SPECIFICATIONS****6.1 DIMENSIONS AND WEIGHT**

Type	MEYRA iMove
Serial number	see nameplate
Weight	approx. 5 kg
Max. load	25 kg
Overall dimensions (L x W x H)	320 x 340 x 90 mm
Frequency	approx. 1 - 12 Hz
Max. amplitude	approx. +/- 4 mm
Travel	approx. 8 mm

6.2 ELECTRICAL DATA

Power supply unit	MW model no. GST60A24-P1J
Input	100 - 240 Volt AC, 50/60 Hz, 1.4A
Output	24 Volt, direct current 2.5 A
Power consumption	max. 60 watt

6.3 MARKING

Certificate	CE
-------------	----

6.4 MATERIAL

The aluminium casing has a robust, impact-resistant powder coating. The remote control consists of anodised aluminium.

7 CARING FOR YOUR iMOVE

Before you start cleaning the iMOVE, disconnect the power supply by pulling out the mains plug. Clean the iMOVE with a slightly dampened cloth and then wipe with a dry cloth. It is essential that no fluid is allowed to enter the iMOVE as this could damage the electronics.

IMPORTANT!

Never use corrosive, abrasive or solvent-based cleaning agents!

8 GUARANTEE POLICY

MotoSana GmbH highly values you as a customer and always endeavours to offer you the best possible service. Every iMOVE device leaves our premises in perfect condition. Should any manufacturing and/or material defects come to light in the course of normal use of the device, we guarantee repair during the first two years from the date of purchase, free of any charge for material and wage costs. The product must have been purchased from an authorised iMOVE partner. Fulfilment of the guarantee assumes transport of the faulty device from the customer to the manufacturer at the customer's expense.

LIMITATIONS IN FULFILLING THE GUARANTEE

The exclusive claim of the customer in the event of a faulty iMOVE is restricted to repair or replacement of the faulty part. Which guarantee service or which combination of guarantee services shall be provided in each case remains at the sole discretion of MotoSana GmbH. Once MotoSana GmbH has established that the product is faulty, MotoSana is entitled to a reasonable period of time for repair or replacement. The replacement product sent to the customer within the scope of the limited guarantee will be manufactured using new and functioning used components. For repaired or replaced products, either the remaining term of the original product guarantee or a guarantee period of ninety days shall apply, depending on which of the two provides a longer remaining guarantee period. The guarantee shall become invalid in the event of a repair being carried out by a person or repair shop not authorised by MotoSana GmbH instead of by MotoSana GmbH itself.

9 PROBLEMS AND HOW TO FIX THEM

- 1 The button on the iMOVE lights up red.
Disconnect the power supply of the iMOVE for five minutes.
If the light remains on, please contact our customer service centre.
- 2 The button on the remote control lights up red.
Replace the battery in the remote control (**required battery: CR2032 button cell**)

SERVICE TELEPHONE
+49 (0)5733 922 - 0

10 MAINTENANCE AND CUSTOMER SERVICE

CUSTOMER SERVICE

Thanks to the use of the best, high-quality industrial goods, no maintenance is necessary. In the unlikely event that a problem does arise, please check the device in line with the possibilities outlined above. Otherwise please contact our customer service centre and have the following information at hand:

- The serial number of your device (to be found on the bottom of the iMOVE)
- Your name
- Your full address
- Your telephone number

IMPORTANT!

Never attempt to open the casing yourself, otherwise the guarantee will lose its validity! No modifications or repairs may be carried out by unauthorised persons.

11 DISPOSING OF YOUR iMOVE

The device must be disposed of in compliance with the relevant national regulations.
There are no environmentally harmful hazardous substances requiring consideration.

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© MEYRA GmbH 2018. All rights reserved.
MotoSana GmbH, Robert-Bosch-Str. 6, D - 86899 Landsberg am Lech, Tel: +49 (0)8191 9406761, info@motosana.de
© MEYRA GmbH 2018. Alle Rechte vorbehalten.
© MEYRA GmbH 2018. Tous droits réservés.
© MEYRA GmbH 2018. Reservados todos los derechos.