



SPECIAL STROLLER

# RACERevo™



CE

ASSEMBLY MANUAL FOR UPHOLSTERY

**AKCES**<sup>®</sup>**MED**



**CAUTION!** Manufacturer takes responsibility for proper usage of the device only when the product was purchased from qualified Akces-Med employee or at specialized medical store.



Medical product first class according to the Regulation 93/42/EWG of European Parliament and Council regarding medical products. For this product the producer issued relevant declaration of conformity. Producer's identification number in the Registration Department of Healing Products, Medical Products and Biocide Products:

**PL/CA01 04088/W**

**Akces-Med Sp. z o.o.** Company introduced quality management system ISO 13485 regarding production, sales and service.

The system is certified by Det Norske Veritas.



# 1. ELEMENTS OF UPHOLSTERY



1. UPHOLSTERY



2.FRAME



3. ABDUCTION BLOCK



4.RAIL



5. CUSHION FOR SEAT



6. HEADREST




7. SIDE COVERS FOR FOLDING MECHANISM



**8. SIDE SUPPORTING PANNELS FOR CALVES**



**9. HARNESS**

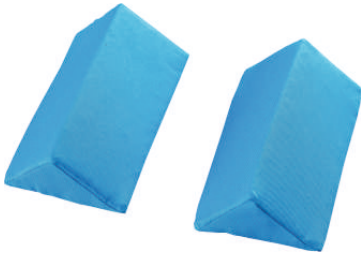


**10. SIDE PELOTTES**



**11. FASTENING BELT**

**12. SPRINGS**

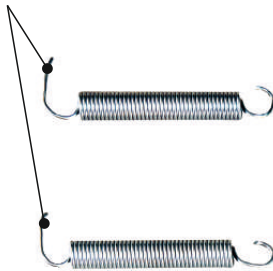

**13. ROD N° 1**



**14. ROD N° 2**



**15. ROD N° 3**



## 2. ASSEMBLY MANUAL FOR UPHOLSTERY

It is necessary to follow these steps so as to mount upholstery on Racer Evo stroller:

**STAGE 1.** Put **UPHOLSTERY (1)** on **FRAME (2)** and backrest and move it backwards.



**STAGE 2.** Insert **ABDUCTION BLOCK (3)** into hole in seat's upholstery.



Screw abduction block under the seat with the aid of knob and washer.



**STAGE 3.** Insert **ROD N° 3 (15)** into hole situated at the front of seat.





Rod should be inserted in a way that its curves are positioned backwards.



**STAGE 4.** Stretch **UPHOLSTERY (1)** and insert curves of **ROD N° 3 (15)** into two holes that are situated on the second side of seat's frame.





**STAGE 5.** Insert **ROD N° 1 (13)** into hole that is situated at the side of stroller. Repeat this activity for the second side.



Rods should be inserted in a way that its curves are positioned backwards.



**STAGE 6.** Stretch **UPHOLSTERY (1)** and insert curves of **ROD Nº 1 (13)** into two holes that are situated on the second side of seat's frame. Repeat this activity for the second side of stroller.



**STAGE 7.** Insert **ROD N° 2 (14)** into hole in upholstery that is situated at the back of upholstery.



**STAGE 8.** Put **SPRINGS (12)** behind **ROD N° 2 (14)** that is situated in hole in upholstery. Take spring and "pierce" upholstery at a distance of 9cm from the end of upholstery on both sides.







**STAGE 9.** If **SPRINGS (12)** are mounted, it is important to bend the endings of springs with the aid of pliers.



**STAGE 10.** Stretch as well as insert the endings of **SPRINGS (12)** into holes that are situated at the back frame of stroller.



**STAGE 11.** Thread **FASTENING BELT (11)** through backrest's frame as it is illustrated in the picture below. The ending of belt should be inserted into regulator. These activities should be repeated for the second ending.









**STAGE 12.** So as to mount **HARNESS (9)** it is necessary to take out four belt that are situated behind upper as well as bottom part of backrest's upholstery.



**STAGE 13.** Thread belts through four clamps in vest.





**STAGE 14.** Place **SIDE SUPPORTING PANNELS FOR CALVES (8)**



**STAGE 15.** Place **HEADREST (6)**



**STAGE 16.** Place **SIDE PELOTTES (10)**



**STAGE 17. Take out CUSHION FOR SEAT (5)**



**STAGE 18. Place RAIL (4)**









ORTHOPAEDIC DEVICES MANUFACTURER



Akces-MED. Sp. z o.o.  
ul. Leszka Czarnego 3  
35-615 Rzeszów  
tel. (17) 864 04 70

[www.akces-med.com](http://www.akces-med.com)

e-mail: [akces-med@akces-med.com](mailto:akces-med@akces-med.com)

Authorised Dealer



Date of issue: 14.11.2014

The last update: 14.11.2014



AKCES-MED Ltd. reserve the right to introduce technical and trading changes  
in the content of the instruction without warning.